



DanyelBIOTECH
at the heart of action



TakaRa




Clontech

Lenti Viral Titration Kits

Lenti-X qRT-PCR Titration Kit

Set up the qRT-PCR reaction with purified RNA from harvested lentiviral supernatant. Determine the viral genome content from a calibrated RNA standard curve.

qPCR

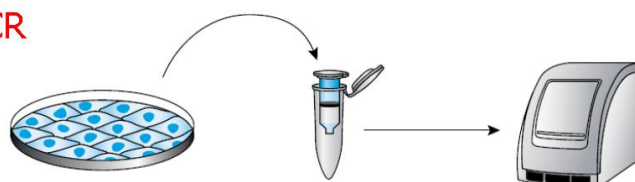

631235


**Lenti-X™ qRT-PCR
Titration Kit**

200 RNX

Lenti-X p24 Rapid Titer Kit

Lentiviral supernatant is bound to anti-p24 coated wells and detected using a combination of biotinylated anti-p24 secondary antibody, streptavidin-HRP, and a development reagent.

ELISA

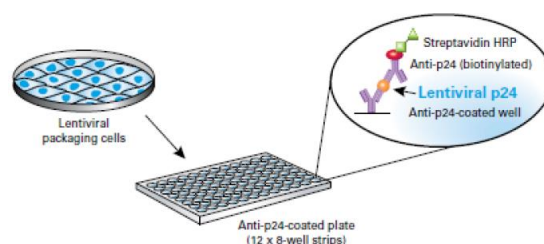

632200

**Lenti-X™ p24 Rapid
Titer Kit**

96 Rnxs

Lenti-X™ GoStix™—Instant Titration

Apply 20 µl of your lentiviral supernatant to the sample well, add 4 drops of buffer, and wait for the indicator bands to appear. Two bands mean you have $\geq 5 \times 10^5$ IFU/ml.

GoStix



631243/4

Lenti-X™ GoStix™

20/50 Rnxs

Please contact us at:
Danyel Biotech: 08-936-6606
Yaniv: 054-4606301
mbc@danyel.co.il

